


# PreciseView System

**Our revolutionary PreciseView precision alignment and positioning system — the most versatile imaging system in the industry — simplifies all your x-ray procedures and reduces your operating costs.**

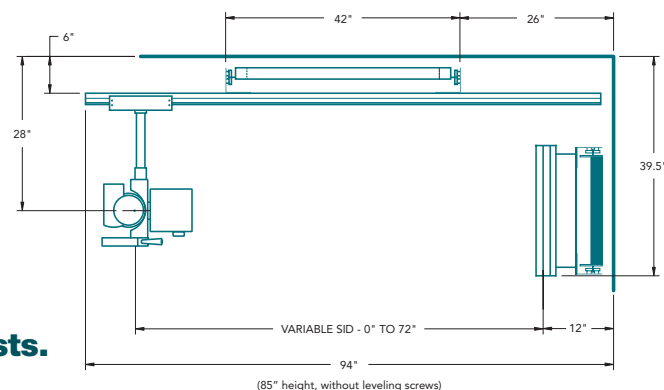

**Universal provides physicians with more generator selections and options than any other manufacturer. Our extensive range of generators are designed and manufactured by Universal, guaranteeing you the quality, reliability and features you demand.**

- **Our generators range** from the budget-friendly 400mA single-phase model to a wide variety of high frequency generators with features like anatomical programming.
- **Universal High Frequency Generators** start at 500mA and provide you with the power to image a full range of patients, dramatically reducing radiation dosage and exposure times resulting in fewer retakes due to involuntary motion.
- **Anatomical programming** allows the operator to select the desired anatomical region and view, then the generator automatically sets the optimal technique factors. Customized programming allows fine tuning the techniques to optimize viewing preferences.

**Our Upright Raymaster System combines economy and flexibility to provide a wide range of diagnostic procedures with a traditional upright system at an affordable price.**

- **Flexible mounting** of the Upright Raymaster tubestand column allows for floor-to-ceiling or floor-to-wall mount interchangeability, for easier, quicker installation.
- **The Upright Raymaster system accommodates** a full range of views and patient heights with smooth, effortless vertical tube travel from 14 1/2" to 69 1/2".
- **Heavy-gauge welded steel** structures provide a sturdy, long-lasting system that will provide many years of reliable use.



**Universal**

[www.delmedical.com](http://www.delmedical.com)

PreciseView Imaging System

**Universal**



**The PreciseView's  
x-ray tube moves on  
a structural steel  
cross-arm to provide  
unparalleled rigidity  
compared to traditional  
floor-to-ceiling/floor-to-  
wall tubestand systems.**

### Ease of operation

While others talk about innovation, for 75 years Universal has demonstrated it. The PreciseView's advanced design delivers a real clinical edge. Aside from offering precise positioning and alignment—a major contributor to high quality radiographs—the PreciseView provides maximum ease of operation with time-saving, single motion positioning for virtually all upright and extremity exams.

### Productivity-enhancing

The PreciseView's engineers relied heavily on input from specialty practitioners to design a system that would streamline operation. It has been designed to make positioning easier for both patient and doctor or technician. Virtually effortless movement is assured by the use of precision stainless steel bearings, which allow smooth and precise tube positioning, both horizontally and vertically. Precise positioning means fewer repeats, less film waste and faster throughput. Single-motion positioning reduces

patient positioning time by up to 80%. This is accomplished through interlocking the tube positioning structure to the image receptor when performing procedures where the central ray is 90 degrees to the patient. Angled views are accommodated by simply disengaging the tube-image receptor interlock through the use of a 'T' handle pin.

e a s e o f o p e r a t i o n  
**image quality**

### Extremity views performed effortlessly

The days of awkward patient positioning for extremity views are over. The PreciseView's unique, tilting image receptor allows a new level of patient comfort in performing extremity views while enhancing the professional image of the practitioner. It's ability to tilt into a horizontal platform and adjust vertically allows positioning equipment to the patient—rather than having the patient conform to the equipment in an awkward and uncomfortable manner. The PreciseView's exclusive and unique features will allow you to accommodate any radiographic study that comes your way.

### A modern system for the modern practice

The progressive look of the PreciseView's counterbalanced structure eliminates unsightly floor and wall rails found on other systems, and makes installation easier. Modular panels contribute to the system's contemporary aesthetic and allow for faster installation and front access serviceability.

### We also thought about the little things

The laser light of the collimator confirms accurate positioning. The collimator exceeds the lead lined protection that you expect for the operator and patient. We also include a built-in spare bulb for quick replacement at those awkward times when needed. Additionally, the collimator is provided with a SID tape measure to facilitate accurate positioning for extremity and angled views.

### A cut above the rest

The PreciseView was engineered for the discerning physician who demands not only high quality films, but also quality equipment that can be used for years to come while enhancing the professionalism of their practice.

### Generators to fit your needs

Choose from a variety of Universal generators—from the basic 400 mA, 125 kVp single phase to state-of-the-art 500 or 600 mA high frequency with or without anatomical programming for use with the PreciseView.

### Five year equipment warranty

Complies with all CDRH performance standards in effect at the date of manufacture.

**Call Universal today for more information and to discuss your digital options.**


**1-800-800-6006**

Made in U.S.A.  
As we strive for constant improvement, specifications are subject to change without notice.  
An ISO13485 Certified Facility.

# Versatile

w i d e r a n g e o f c a p a b i l i t i e s

**The PreciseView's fully counterbalanced wallstand is supplied with a 17" x 17" grid cabinet and 103 line, 8:1 ratio grid. Extensive travel will accommodate a full range of views and patient heights.**

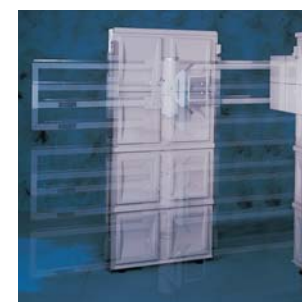

# Reliability

a c c u r a t e d i a g n o s i s